

**AMIRO**

**AMIRO U1 UltraLift Skin-lifting Ultrasound Home Device**



**User Manual**

# AMIRO

\* We appreciate you buying an AMIRO product.

- For your safety, please read this manual carefully before use, and strictly follow the instructions provided.
  - Please save this manual for future reference.
  - Please include this manual when transferring ownership of the product to others.

## **Instructions for Use**

Dear customer, thank you for choosing AMIRO U1 UltraLift Skin-lifting Ultrasound Home Device from Shenzhen Accompany Technology Co., Ltd. (“our company”).

Before using the device, please ensure that you have carefully read and fully understood the user manual. To use the device properly, make sure the instructions are followed carefully.

If you have any questions while using the device, feel free to reach out to our company for assistance.

We promise to offer you the best of our service at all times.

## **Safety Warnings**

- When using the product, please strictly follow the instructions in the manual.
- Device adjustment and maintenance should be carried out by designated and trained personnel of our company.
- The electrical equipment and work environment must meet national standards, industry norms, and the manual requirements.
  - Strictly adhere to the instructions for the use, charging, and storage of the device.
  - Use only the accessories provided by our company

## Contents

<b>Device and Manufacturer Information</b> .....	2	<b>Main Technical Parameters</b> .....	29
<b>Safety Instructions</b> .....	3	<b>Function description</b> .....	30
<b>Contraindications and applicable scenarios</b> .....	13	<b>Product maintenance</b> .....	31
<b>Operating principle</b> .....	15	<b>Product storage and transport environment</b> .....	31
<b>Product Components</b> .....	17	<b>Product troubleshooting</b> .....	32
<b>Product Features</b> .....	18	<b>Main safety properties</b> .....	34
<b>Product Usage</b> .....	19	<b>Explanation of graph and symbol</b> .....	34
<b>How to Use the Product</b> .....	21	<b>Waste disposal declaration</b> .....	35
<b>Product Instructions</b> .....	22	<b>Hazardous Substance Label</b> .....	36
<b>Mode Description and Switching</b> .....	24	<b>Warranty Policy and Warranty Certificate</b> .....	37
<b>After Skincare</b> .....	28		

## **Device and Manufacturer Information**

---

### **Product Information**

Name: AMIRO U1 UltraLift Skin-lifting Ultrasound Home Device

Model: ABU001

### **Manufacturing Information**

Manufacturer: Shenzhen Accompany Technology Co., Ltd.

Address: F26, High-Tech Zone Union Tower, 63 South 10th Avenue, High-Tech Zone, Binhai Community, Yuehai sub-District, Nanshan District, Shenzhen, China

Postal Code: 518063

Contact number: 400-048-0888

After-sales Service Provider: Shenzhen Accompany Technology Co., Ltd.

Address of the After-sales Service Provider: F26, High-Tech Zone Union Tower, 63 South 10th Avenue, High-Tech Zone, Binhai Community, Yuehai sub-District, Nanshan District, Shenzhen, China

## Safety Instructions

---

To use this device safely and to prevent bodily harm to the user or others, please read the safety precautions carefully and follow the precautions outlined.

About the symbols and their meanings:



Warning: Operations that could potentially harm the user or damage the device



Caution: Additional instructions are provided for better device usage.



### **Warning: Use, Charging and Storage Safety**

To prevent injury or discomfort, and to ensure the device's life span, please ensure that the environment where the product is used and stored is dry, free from flammable or explosive materials, not exposed to strong direct light, and kept at an appropriate indoor temperature.

- Keep away from fire sources or flammable materials during use, charging, or storage.
- Avoid proximity to strong electromagnetic fields, shortwaves, or microwaves to prevent interference and unstable energy output.
- Refrain from using, charging, or storing device in damp or immersion-prone areas to prevent electric leakage or malfunction.
- Avoid using, charging, or storing device in environments with high levels of corrosive substances (such as hydrogen sulfide, sulfur dioxide, ammonia, nitrogen, hydrochloric acid, sulfuric acid, etc.) to prevent corrosion of the main unit and accessories
- Avoid using, charging, or storing the device in areas with significant temperature fluctuations (near heaters, hot pots, air conditioning, etc.) to prevent potential risks due to varying temperature.
- To prevent dropping or damaging the device, avoid using, charging, or storing it in tilted, unstable, or vibrating conditions.
- To protect the environment and prevent skin or body discomfort, please prevent the device from coming into contact with pollutants or waste during use.



### Warning Biological Hazards and Contact Safety

- The components of this device that come into contact with the human body are crafted from high-purity aluminum alloy, non-toxic PCTG materials, and 17-4 stainless steel. These materials have successfully undergone nickel release (GB/T 19719-2005 & GB/T 28485-2012) and biocompatibility testing (GB/T 16886.5-2017 & GB/T 16886.10-2017), ensuring they are safe for use. If you have a history of allergies to these materials, please avoid direct contact with the device or stop using it if you experience any allergic reactions. Consult a doctor if necessary.
- This device is for personal use on the face only. Please refrain from sharing it with others. To avoid cross-infection when the device is shared among multiple users, remember to clean the surface of the ultrasonic energy head before each use.
- Due to the presence of bacteria on facial skin, it is important to clean the electrode head thoroughly before and after each use to prevent bacterial growth.
- For hygiene purposes and to avoid cross-infection of bacteria or viruses, do not use this device on non-human animals, such as pets.




### Warning: Hazards of Water Damage

Classified as a Class III electrical device, the device has successfully passed IPX6 waterproof testing. It must be connected to a power source for operation and charging.

- Prevent liquid from entering the main unit or charging port. After cleaning, ensure the device is completely dry before storing it in a dry environment to avoid potential short circuits, electrocution, or fires caused by water damage.
- To prevent electrocution, short circuits, or fires, refrain from using, charging, or storing the device near water sources, in bathrooms, or while bathing.
- Do not immerse any part of the device in liquid to prevent electrocution, short circuits, or fires.
- The device's bottom charging port features a waterproof plug design. When cleaning, make sure the waterproof plug is securely in place to prevent liquid from entering the charging port and potentially causing electrocution, short circuits, or fires.



✓ Waterproof plug installed correctly.

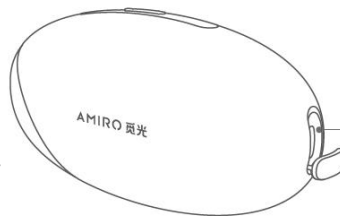

✗ Waterproof plug not installed.



**Warning: Stamping (intense mode, gentle mode) should not be applied to the following areas:**

- Avoid using this device near the eye socket and eyeball to prevent tissue damage.
- To prevent tissue damage due to depth concerns, do not apply this device within 1.5 cm in front of the tragus, in the philtrum, around the mouth corners, in the nasal wings, within the throat, along the jawline (mandible), and within the thyroid and laryngeal regions.
- Please note that this device is intended to only be used on the face and should not be used on other parts of the body.
- Use of the device with a facial mask in these areas is not recommended for users with a small face.

Avoid using this device on the area around the eyeball as well as upper and lower eyelids.

Do not apply this device on the philtrum area.

Avoid using this device on the trachea and lymphatic area of the neck.

Maintain a distance of at least 1.5 cm in front of the tragus when using the product.

Avoid stamping the same area repeatedly and make sure to follow the grid sequence for stamping.

Do not use this device on the edge of the lower jaw (mandible) area.



**Warning: Precautions for the use of different stamping modes (intense mode, gentle mode)**


**Face Mask**

- When using this device for stamping treatment (Intense/Gentle Mode), it should be paired with AMIRO U1 UltraLift Timeless Four-zone Face Mask. Avoid using other face masks to prevent potential injury.
- The face mask is designed for single use only. Avoid recycling, washing, drying, or attempting to reuse it.
- When the stamping area of the face mask appears dry, white, or lifted, replenish the essence/gel, thoroughly moisten it, and proceed with the treatment.




## Face mask + device

- When using this device for stamping treatment (Intense/Gentle mode), ensure that the device head's surface fully contacts both the face mask and facial skin. Avoid partial contact with the skin. [Figure 1]
- When stamping, make sure to follow the grid pattern precisely. [Figure 2, Figure 3]
- The device should be positioned with the smart sensor facing upwards when using the vertical grid. [Figure 4, Figure 5]
- The device emits energy once it detects the color blocks on the mask. To prevent potential injury, stamp the device steadily without moving or shaking it.
- Avoid continuous stamping. Adhere strictly to the grid sequence for stamping.
- Follow the recommended frequency of use: It is recommended not to use stamping treatments more than two times a week (including: Intense & gentle mode).

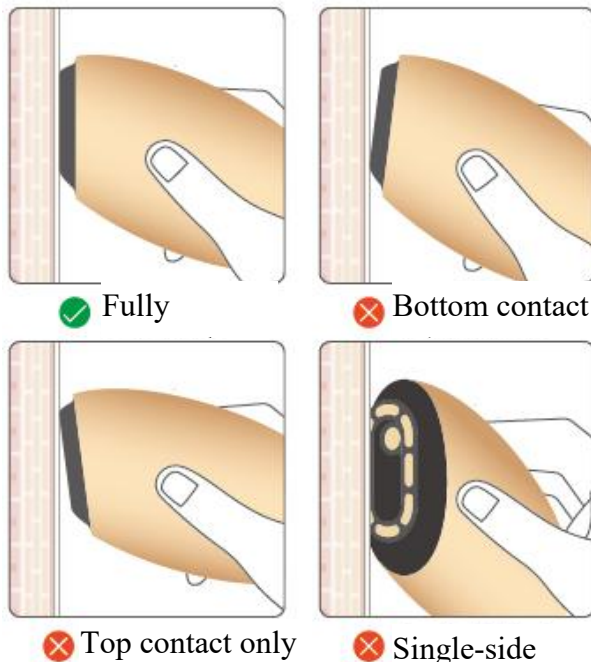

[Figure 1]

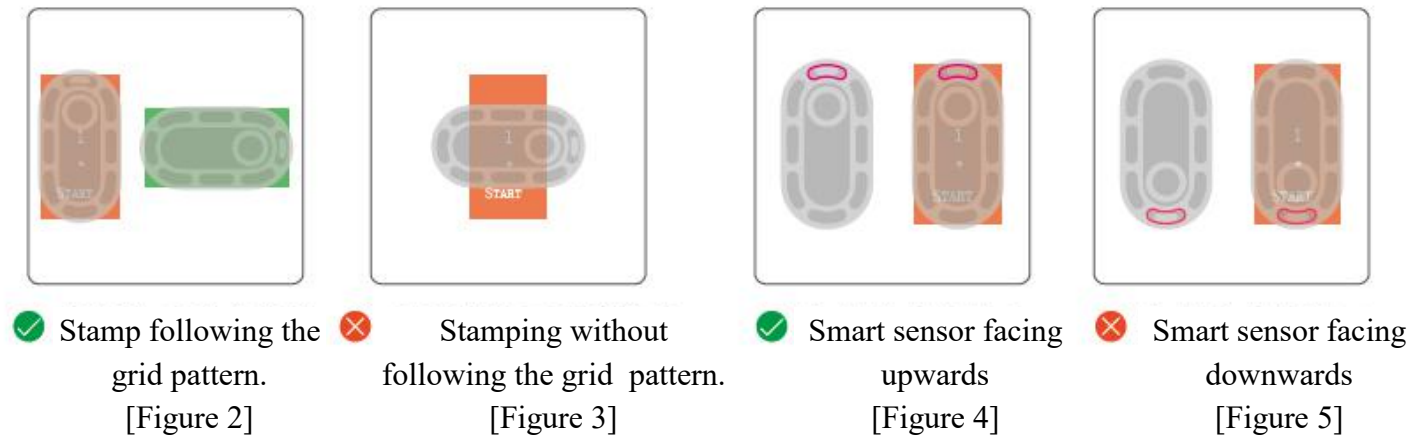



### Warning: Areas to be avoided for Express Mode

- Avoid using this device near the eye socket and eyeball to prevent tissue damage.
- To prevent tissue damage due to depth concerns, do not apply this device within 1.5 cm in front of the tragus, in the philtrum, around the mouth corners, in the nasal wings, within the throat, along the jawline (mandible), and within the thyroid and laryngeal regions.
- Please note that this device is intended to only be used on the face and should not be used on other parts of the body.




Avoid using this device on the area around the eyeball as well as upper and lower eyelids.



Do not apply this device on the philtrum area.




Avoid using this device on the trachea and lymphatic area of the neck.

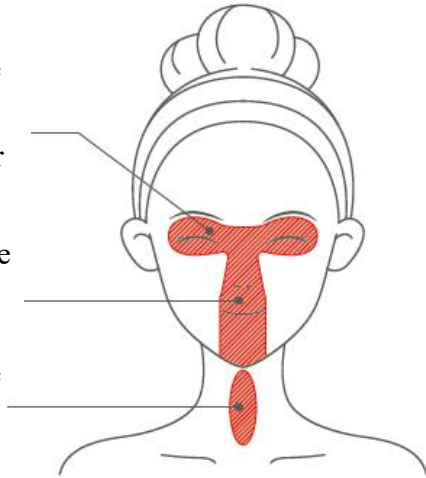

Maintain a distance of at least 1.5 cm in front of the tragus when using the product.



Avoid gliding or circling the device in the same area repeatedly

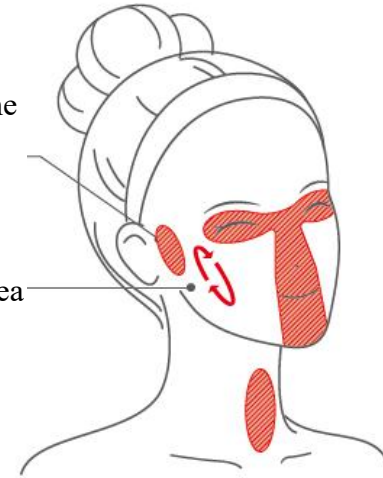



**Warning: Instructions for using the gliding mode (Express mode)**

- When using this device for gliding treatment (Express mode), ensure consistent contact between the device head and the face. Avoid any abrupt movements of the head during contact to prevent potential injury. [Figure 1]



✓ Gliding over the face with appropriate amount of contact.



✗ Do not shake the device head after contact.

[Figure 1]



**Warning: Risks of improper operation**


To ensure the effectiveness and safety of the product, make sure that all functions and components of this product are in good condition before use.

- Do not use the device for non-cosmetic purposes.
- Do not use it for pet care to prevent injuries to your pet.
- To prevent physical discomfort and skin damage, avoid overusing it or similar devices.
- Remember to switch off the power after using the device or when not in use.
- Use the device after you have removed your makeup.

- Feeling different temperatures while using the device is normal, as individual skin sensitivity differs. Adjust the temperature to your preference as necessary.
- To prevent skin damage or discomfort, it is recommended to avoid prolonged use on the same area.
- The ultrasonic energy head contains a liquid that is essential for the stable conduction of ultrasonic waves. Do not remove the film covering the energy head surface. Inspect the head for any signs of damage or leakage before use. [Figure 1]
- To avoid irreversible deformations or alterations to the film covering the energy head's surface, which could negatively impact ultrasonic energy output, avoid touching the film with sharp objects. [Figure 2]



[Figure 1]

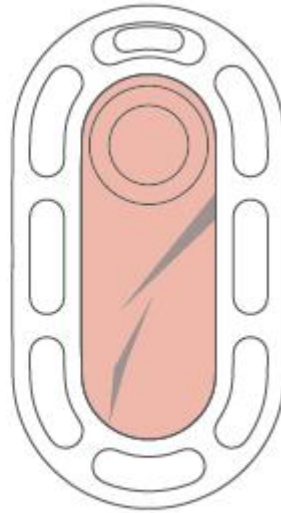

[Figure 2]

Do not use the product for purposes other than those listed in the manual and avoid the following operations to prevent personal injury or device failure.

- Avoid using the device near strong magnetic or high-frequency devices such as microwave ovens.
- Avoid contact with sharp objects that may lead to device malfunction.
- Avoid applying the device on injured or healing skin; refrain from using the device when skin conditions are present, whether on the entire body or on specific areas.
- To prevent device malfunctions and personal injury, ensure that the electrode pads do not touch metal while in use or during charging.

- To prevent component failures or injuries from dropping the device, make sure to keep the main unit and base secure and stable to avoid any impact or falls.
- Do not use the device if it is malfunctioning as doing so may result in bodily injury or danger.
- To avoid overheating hazards, do not cover the main unit with materials that can hinder heat dissipation, like fabric or blankets.
- Avoid putting excessive weight on the main unit or pulling on the charging cable.
- Keep the device away from children as it contains small parts that could be swallowed or cause choking.
- Properly dispose of plastic packaging bags after unpacking to prevent risks of children suffocating or pets playing with them.
- When using, charging, or storing the device, make sure to keep it out of reach of infants, toddlers, and individuals with mental disorders to prevent any potential injuries.



### **Warning: Electrical Safety**

- Do not disassemble, repair, modify, or damage the main unit or charging cable. Do not bend, pull, or twist the power cord.
- Do not expose the main unit or charging cable to water or high humidity to prevent electrocution or leakage. Avoid charging, using, or storing the device in areas prone to moisture, such as bathrooms or sinks.
- This device does not come with an adapter. When charging, make sure to use an adapter that complies with national certification standards and provides a 5 V = 3 A output.
- Do not connect the device to unstable power sources to prevent inconsistent charging efficiency or damage to the device.
- Do not connect this device in parallel or series with other electrical appliances. When the rated capacity is exceeded due to improper usage such as connecting multiple devices in parallel, it may lead to a fire hazard.
- Only use the charging cable provided by our company. Using other cables may cause damage to the main unit.
- Avoid using improper adapters to prevent safety risks like overheating, electrocution, and leakage.



### **Warning: Risk of Electrocution**

- Do not use the device if you are wearing contact lenses or metal accessories such as glasses, earrings, necklaces, watches, or rings.
- To avoid risking electrocution, do not attempt to open the main unit's casing by yourself.
- Do not handle the charger with wet hands. Ensure that the device, charging cable, and adapter are dry at all times.
- Do not use the charging cable if it is loose or damaged when inserting it into the adapter.
- Hold it properly when pulling it out. To avoid the risk of electrocution, do not plug or unplug the charging cable with wet hands. After removing the charging cable, ensure that the waterproof plug is secured in place.
- Do not use unapproved adapters, as they may pose safety risks such as overheating, electric shock, or power leakage.
- If the charging cable is damaged, contact the manufacturer or dealer for a replacement.



**Caution: Improper storage may reduce the lifespan of the device**

Proper storage is crucial for the effective use of this product. Position it on a secure surface and ensure the safeguarding of the device head.

- Keep the main unit and charging cable out of reach of children, infants, or pets.
- Do not place the device in a microwave or high-pressure container.
- Disconnect the charging cable after using the device or while cleaning the main unit, during power outages, lightning, or similar situations to avoid any potential damage.
- Prior to use, inspect the device head for any cracks, deformations, or damage. Avoid use if any abnormalities are found.
- Do not submerge the device in liquid-filled containers.
- Do not clean the device with solvents like alcohol, diluents, gasoline, nail polish remover, diesel, makeup remover, etc.



**Caution: Transportation Safety**

- When transporting the device, avoid shaking or placing heavy objects on top of it. Otherwise, it may lead to damage or malfunction of the device.
- When transporting the device, use proper packaging and handle it carefully to prevent the main unit from falling or experiencing heavy impacts that may cause damage.
- Avoid exposing the device to rain or submerging it into water during transportation to prevent the risk of short circuits or electrocution.

## Contraindications and applicable scenarios

---

People with any of the following conditions may not use this device as it could result in accidents or skin and bodily injury:

- Experiencing allergic reactions.
- Skin areas with impaired integrity.
- Skin areas with damage, acute inflammation, active infections, rosacea, ongoing vitiligo, or severe acne.
- Suffering from severe systemic illnesses.
- Areas with metal implants or fillers.
- Undergoing wound treatment or managing systemic or local skin conditions.
- Please consult a doctor before using the device if you have had permanent dermal fillers injected or have received temporary dermal fillers/Botox injections in the facial areas of application.
- Individuals with abnormal blood pressure, heart disease, cardiac functional disorders, or tendencies towards functional disorders.
- Infectious diseases or fever.
- Hemorrhagic disease or cicatricial diathesis.
- If you have thyroid nodules or other thyroid disorders, consult your physician before using this product.
- Individuals with abnormal pain perception, inability to express consciousness, or lack of temperature and stimuli perception.
- If you are undergoing dental or oral treatment, consult a doctor before using the device.
- Individuals experiencing menstruation (which can lead to hormonal instability and significant skin changes), pregnancy, or breastfeeding.
- Infants, young children, or individuals under 18 years old (in adolescence).
- Individuals using implantable electronic devices such as pacemakers.




- Individuals using life-sustaining electronic devices such as heart-lung machines.
- Individuals using wearable medical devices such as electrocardiographs, infusion pumps.
- Suffering from, or having previously suffered from, cancer (especially skin cancer and precancerous moles), acute illnesses, tuberculosis, febrile diseases, infectious diseases, blood diseases, visceral diseases, skin diseases, or allergic dermatitis
- Impaired immune system due to immunosuppressive diseases such as AIDS or the consumption of immunosuppressive drugs
- Experiencing fatigue, illness, under the influence of alcohol or on sleeping pills
- Diseases caused by thermal irritation within the application area, such as recurrent herpes simplex
- Consult a doctor before use if you have high myopia (above 600 degrees).
- The device could cause damage to the skin if it is used repeatedly or for a long period of time on the same area.
- Avoid bypassing or removing the built-in safety mechanism.

Safety Instructions: The list of safety measures above is not exhaustive. Should you have any questions regarding the use of this product, please do not hesitate to contact us or consult a physician.

### **Applicable scenarios**

Target Consumers: Intended for individuals aged 18 and above with the following facial concerns: (Please refer to the information regarding contraindications)

- Undefined facial contours
- Weak jawline
- Loose facial skin
- Target areas: Facial region (please refer to  Warning: for information on the prohibited areas for different modes)

### **Expected usage environment**

Expected usage environment: Suitable for home use.

## Operating principle

---

### Ultrasound

Using high-energy ultrasonic technology, the Amiro AMIRO U1 UltraLift **penetrates the epidermis and targets the SMAS layer**, causing skin, collagen cells, and collagen fibers to contract. Through this mechanism, new collagen fibers are produced and the face is reshaped, resulting in a face lifting effect.

### Micro-current


The principle behind micro-current is that pulse current can stimulate muscles to contract, helping with facial lifting and skin firming.

### Ingredients:

- **When combined with the AMIRO U1 UltraLift Timeless Four-zone Face Mask**, collagen, fibronectin, elastin, and hyaluronic acid sodium are introduced, which are vital components for replicating the dermis layer composition of the human body. With this blend, your skin is tightened, wrinkles are diminished, and phospholipids are incorporated to encapsulate highly active caffeine to ensure precise penetration, thereby improving your skin's microcirculation. Additionally, L-Ergothioneine curbs the expression of TNF- $\alpha$  and MMP-1 in dermal fibroblasts, combats free radicals, and decreases lipid peroxidation.
- **In conjunction with the AMIRO U1 UltraLift Timeless Face Lifting Gel**, low molecular weight collagen penetrates deeply into the skin, stimulating collagen production. As L-ergothioneine enters cells and mitochondria via OCTN-1 transport, it eliminates reactive oxygen radicals and provides antioxidant protection for cells. Furthermore,  $\beta$ -glucan, which has a low molecular weight, is easily absorbed through the skin, refines the texture of the skin, and stimulates the development of new collagen, while CHONDRUS CRISPUS extract enhances keratinocyte differentiation, increasing the moisture content of the skin.


## Enhancing penetration: 1+1>2

Through replicating the human dermis layer composition and providing essential elements such as collagen, fibronectin, elastin, hyaluronic acid sodium, and more to the skin, this treatment stimulates collagen synthesis synergistically, tightening the skin and refining facial contours as a result.

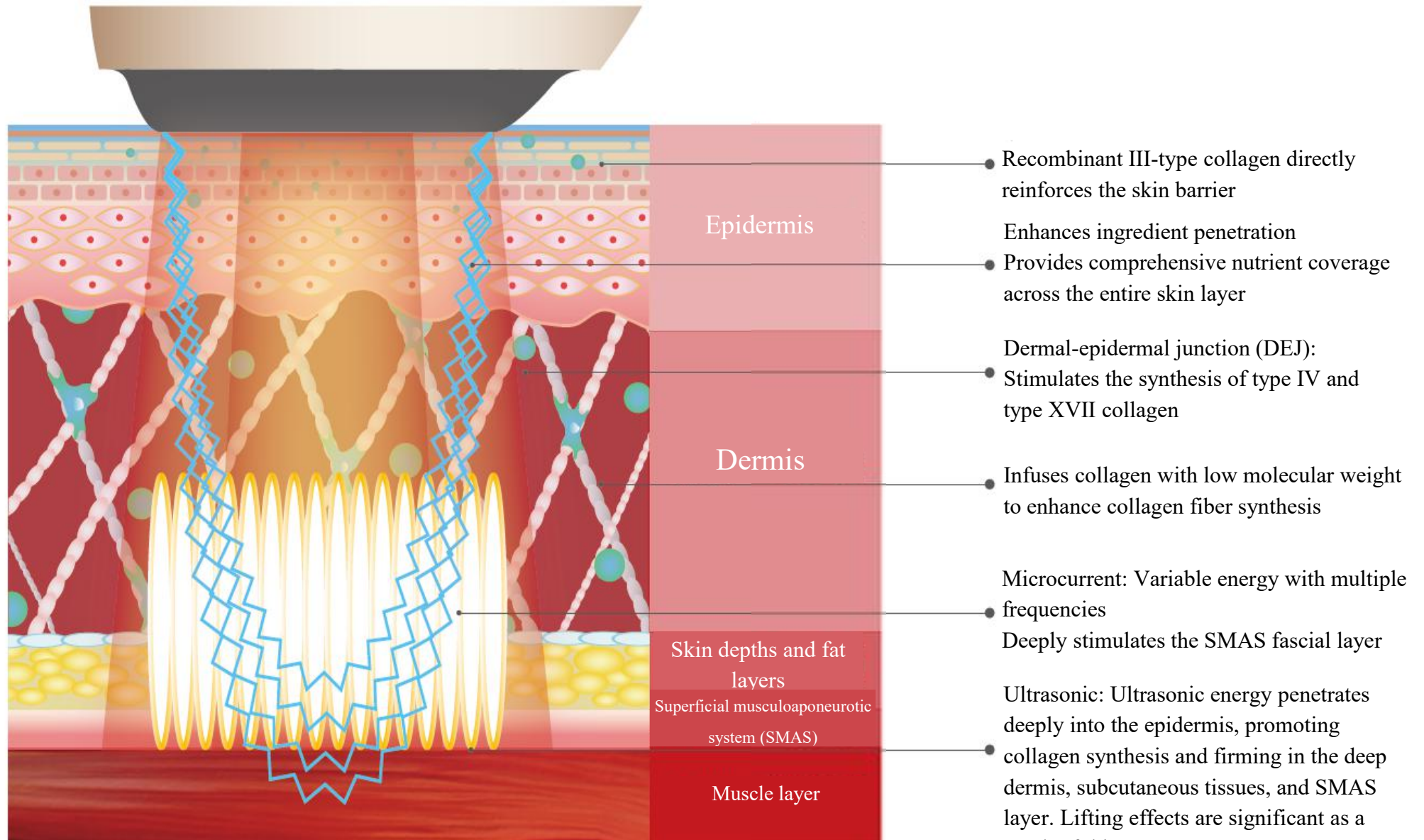


## Product Components

### Product Components

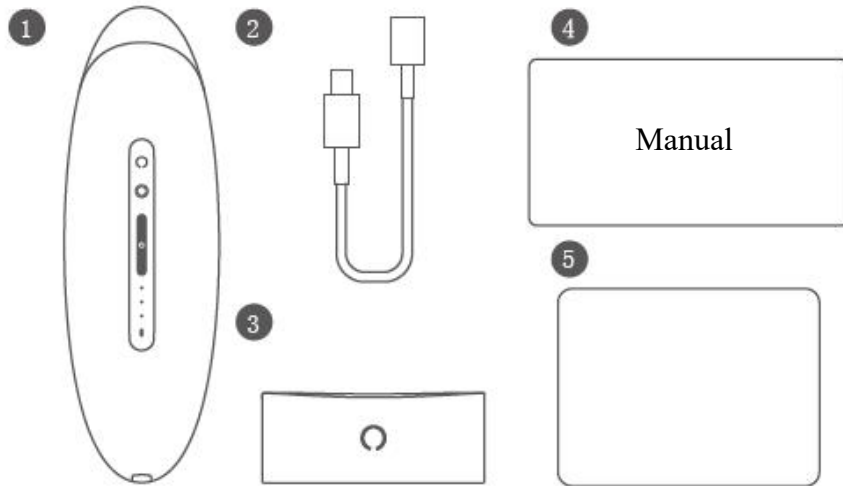


- 1 AMIRO U1 UltraLift Skin-Lifting Ultrasound Home Device \*1
- 2 Charging Cable ×1
- 3 Charging Base\*1
- 4 User Manual & Instructions
- 5 Cleaning Pad\*1

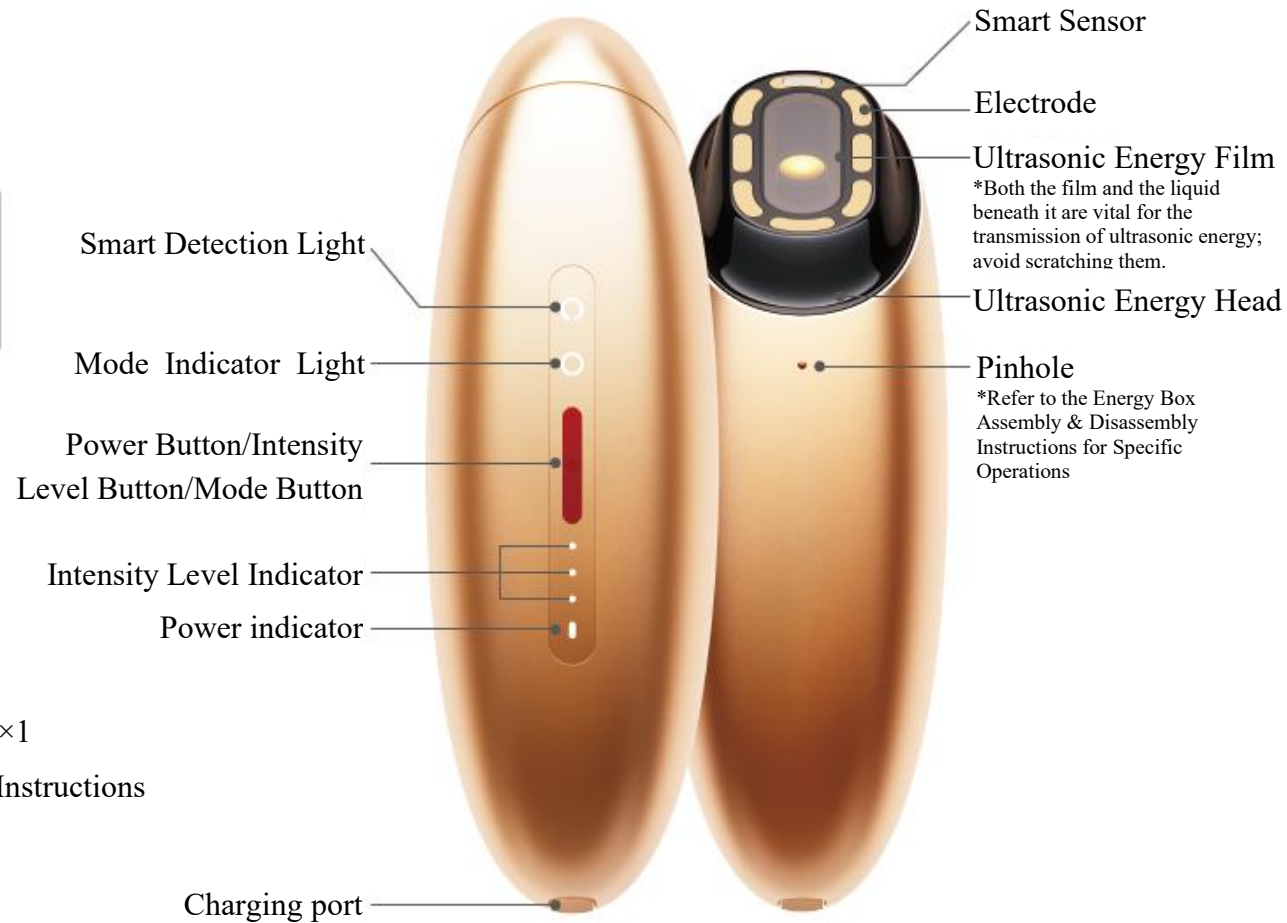





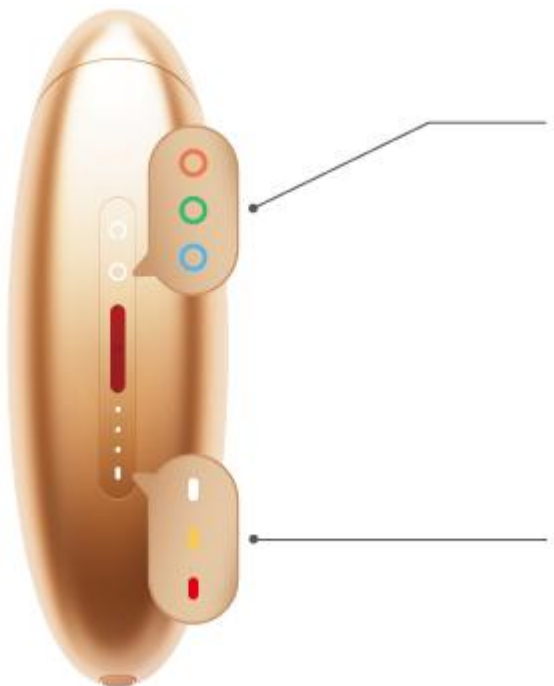


## Product Features

### Power indicator

When the device is turned off, press the red power button briefly to check the battery level



Mode/Function Description		Amiro Ultrasonic Facial Skin-Tightening Machine
○ Intense Mode	3 intensity levels:	Ultrasonic + Microcurrent + Intelligent Zone Recognition + Penetration Promotion
○ Gentle Mode	3 intensity levels:	Ultrasonic + Microcurrent + Intelligent Zone Recognition + Penetration Promotion
○ Express Mode	3 intensity levels:	Ultrasonic + Microcurrent + Penetration Promotion

Battery Level Indicator	Charging Status Indicator
 If the red light remains on for an extended period, the battery is low and insufficient for one use; please charge the device.	 Red light flashing: Charging, critical battery level
 Yellow light on: Low battery level, please charge after use	 Yellow light flashing: Charging, low battery level
 White light on: Battery partly charged	 White light flashing: Charging, battery partly charged   Green light on: Battery fully charged, unplug the charging cable and secure the waterproof charging plug cover.

## Product Usage

---

### Prep and Safety Tips for First Use

Open-box checking: Unpack the device and its accessories, cross-check them against the packaging list, and keep the packaging materials for transport, storage, and reference in the future.


Check the device and accessories for any damage. If there are signs of mechanical damage or equipment malfunction, do not use the device and contact the distributor or after-sales service immediately.

### Charging

Before using the device for the first time, charge it until the battery indicator shows a steady green light, indicating that the charging process is complete.

#### **Caution:**


- How should the energy head surface come into contact with your skin

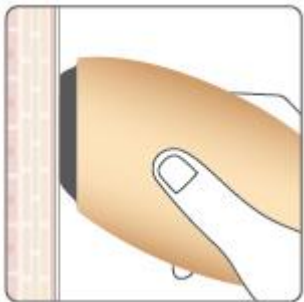

 Proper full contact



 Bottom contact only

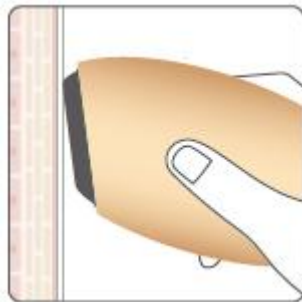



 Top contact only



 Single-side contact only

**⚠ Caution:**

- Do not use the device continuously in the same area.
- If a strong tingling sensation is felt on the skin, lower the intensity or move to a different area.
- A mild tingling sensation near the cheekbones, jawline, or teeth during use is considered normal and safe.
- Use the specific skincare products supplied with this device.
- Avoid using the device on the indicated facial areas highlighted in red:

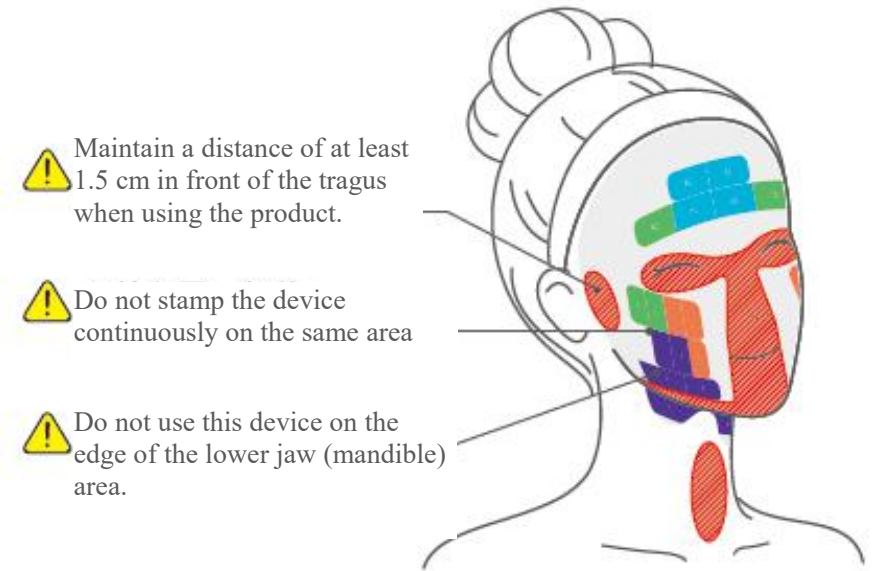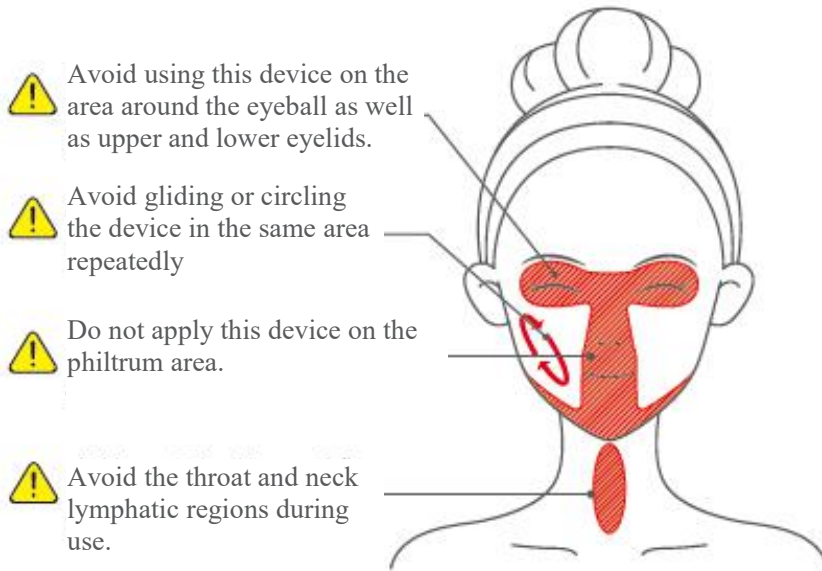


## How to Use the Product




### 1. Cleanse your face

After washing your face, ensure it is thoroughly dry. Remove any jewelry worn on the face, ears, neck, and hands, as well as contact lenses.

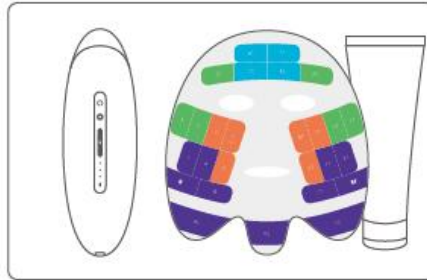

### 2. Prepare the product

Press the button once to check the battery level. If the battery indicator is red, please charge the device before use.

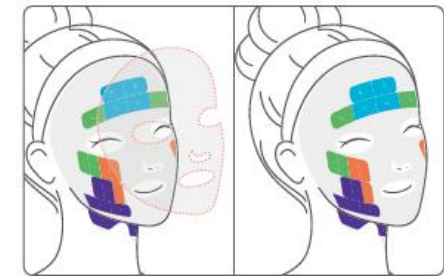


### 3. Apply the face mask

Take the mask out of its package. Place the mask on your face with the **lining paper facing outward**. Peel off the lining paper and delicately press the mask onto your face.



### 4. Switch on the Device

To turn on the device, press and hold the red power button until you feel a vibration. If this is your first time using the device, we recommend that you start at Level 1 and adjust the intensity gradually by clicking the button.

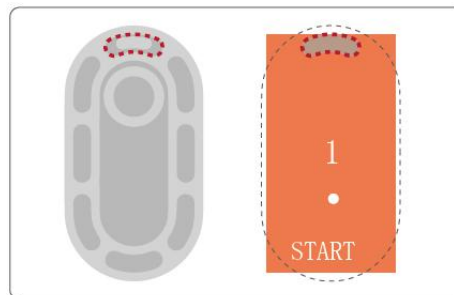

### 5. Start stamping

Place the **smart sensor of the device head securely in the designated slot**. Each successful stamp will trigger a vibration indicating a zone change.



### 6. Stamp sequentially

Stamp three times in a specific order, avoiding repeating stamping on the same area more than three times in a row. The mask should be worn for 10-15 minutes following stamping.



## Product Instructions

---

### Switch on Device


Press and hold the red power button

### Mode Switching

Double-click the red power button


### Intensity Level Switching

Click the red power button once

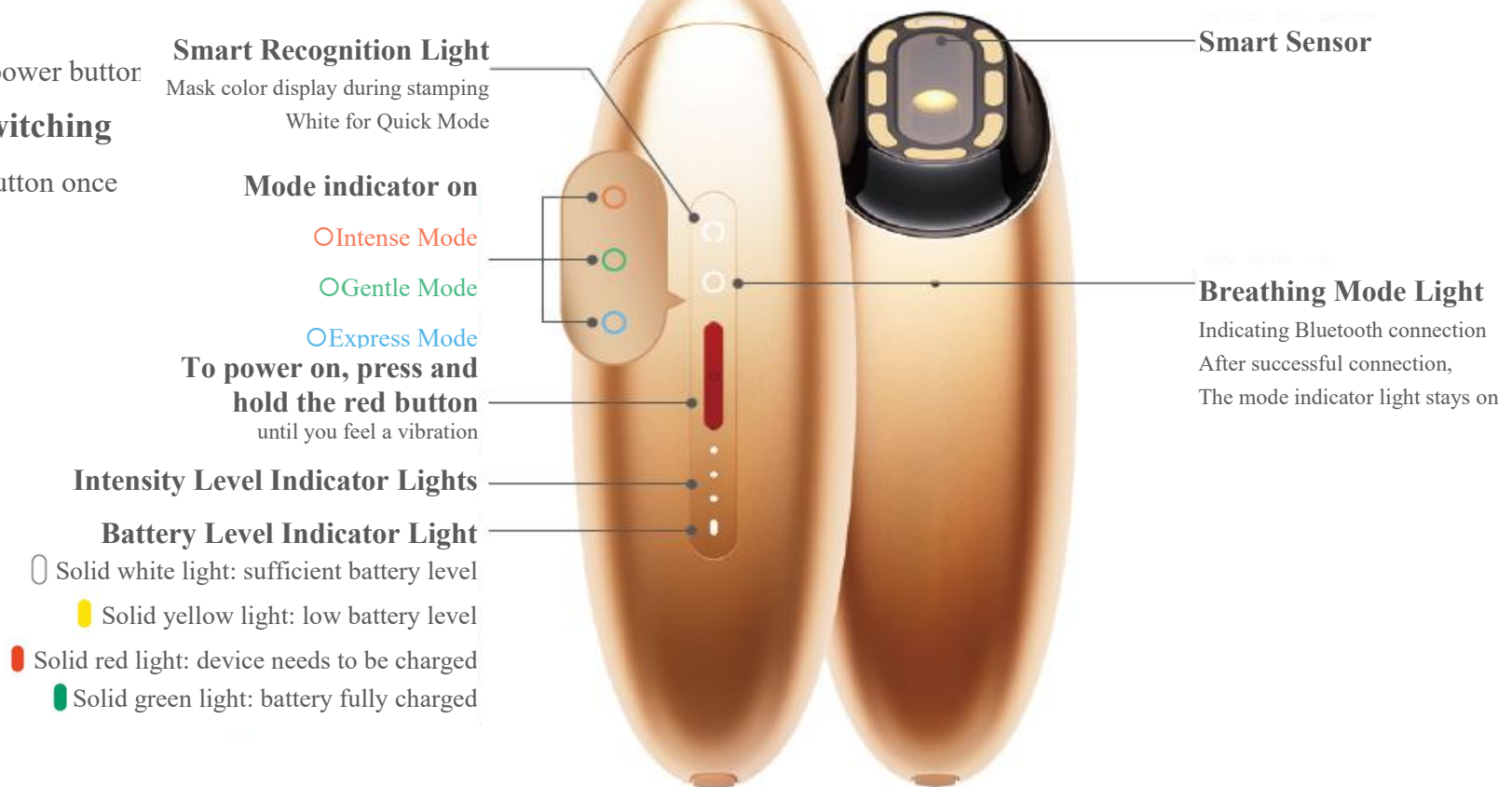

## Select Mode

Double-click to switch between modes: **Intense Mode** → **Gentle Mode** → **Express Mode**

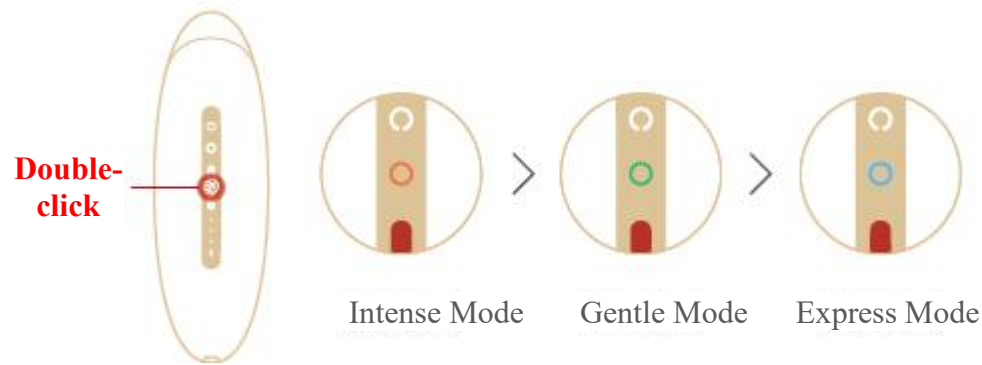


\* Once treatment for a certain mode is finished, the device will vibrate and shut down. To select a different mode, switch on the device again and choose your preferred setting.

## Select Intensity Level



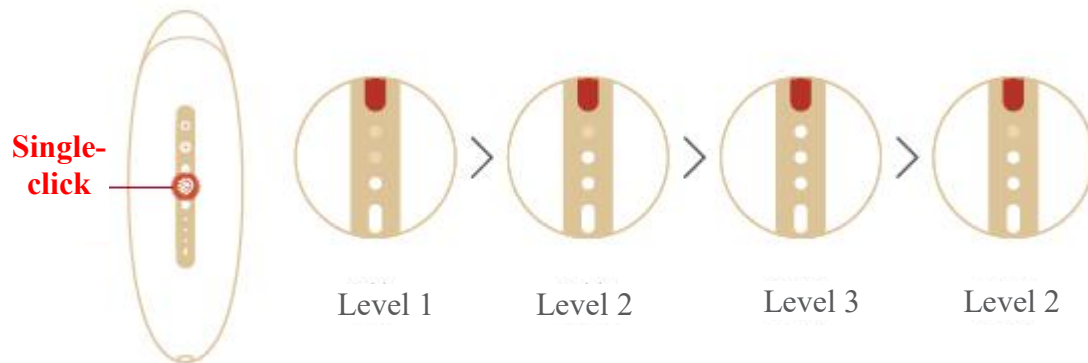

Single-click to switch between intensity levels: **Level 1**→**Level 2**→**Level 3**→**Level 2**→**Level 1**


**Smart Sensor:**



## Mode Description and Switching

### [Intense Mode]

<b>Matching skincare product</b>	Timeless Face Lifting Mask		
<b>Mode Indicator Light:</b>	Orange 		
<b>Smart Recognition Light</b>	Forehead		Eye Contour
	Middle area of the face		Lower Jaw
			
<b>Treatment Areas:</b>	Whole Face		
<b>Matching Methods</b>	Stamp sequentially according to the numbers. Stamp each section three times.		

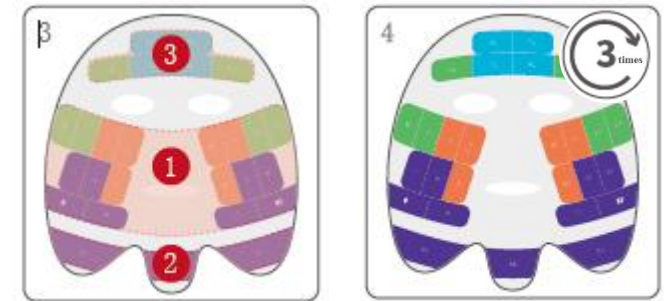

Skincare routine for the face, jawline, and forehead in order


Stamp the whole face for 3 times



Intense Mode: orange light



Timeless Face Lifting Mask

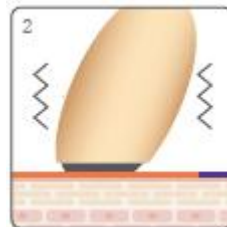

After 5 seconds, the device will vibrate to signal advancement to the next grid.

### Usage Frequency

1-2 times per week


### Warnings

1. The essence features a dense texture to ensure optimal energy transmission with minimal loss
2. Due to the characteristics of ultrasound technology, the Intense Mode requires the official mask that matches the device.
3. **By double-clicking, you can switch to Gentle Mode**, which is recommended for unstable skin conditions.
4. Do not stamp continuously on the same area.

[Gentle Mode]


**Matching skincare product** Timeless Face Lifting Mask

**Mode Indicator Light:** Green 

**Smart Recognition Light** Forehead  Eye Contour   
Middle area of the Lower Jaw   
face 

**Treatment Areas:** Whole Face

**Matching Methods** Stamp sequentially according to the numbers. Stamp each section three times.

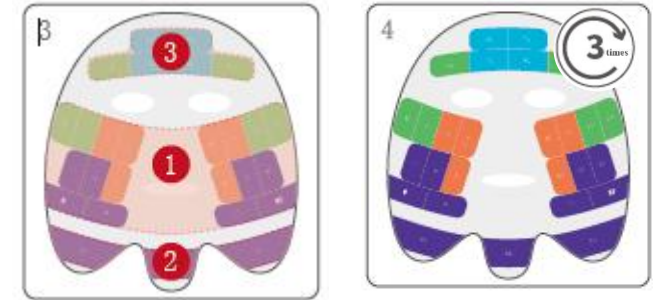

Skincare routine for the face, jawline, and forehead in order


Stamp the whole face for 3 times



Gentle Mode: green light



Timeless Face Lifting Mask

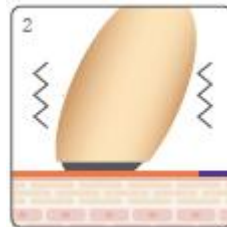

The device vibrates for 5 seconds to indicate advancement to the next grid.

**Usage Frequency**


1-2 times per week

**Warnings**



1. The essence features a dense texture to ensure optimal energy transmission with minimal loss
2. Due to the characteristics of ultrasound technology, the Gentle Mode requires the official mask that matches the device..
3. Avoid repeatedly stamping on the same area continuously.

## AMIRO BEAUTY U1 UltraLift Timeless Four-zone Face Mask
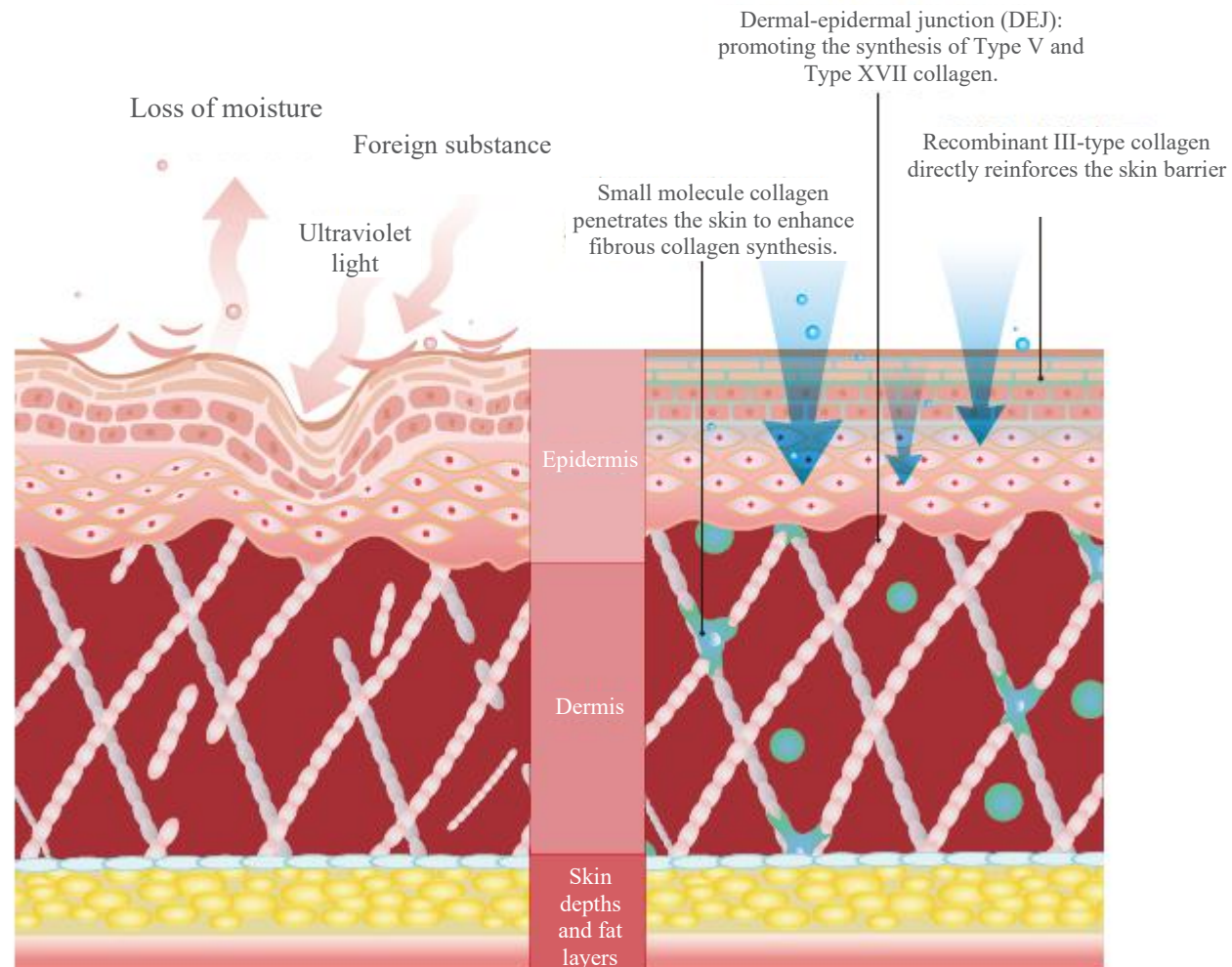
### Working principle illustration:




[Express Mode]


<b>Matching skincare product</b>	Timeless Face Lifting Gel
<b>Mode Indicator Light:</b>	Blue 
<b>Smart Recognition Light</b>	White 
<b>Treatment Areas:</b>	Whole Face
<b>Matching Methods</b>	Glide over the entire face. Vibration of the device indicates zone change



Quick Mode: blue light

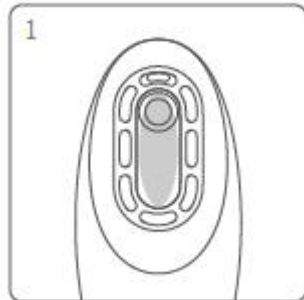


Usage of gel/solution:

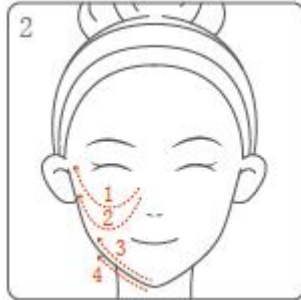


**Left Cheek: 50 seconds**

- Apply gel/solution on the left cheek
- Repeat Steps 1-4 (1 set) 3-4 times
- After 50 seconds of use, the device will vibrate briefly, prompting a zone change to the right cheek.

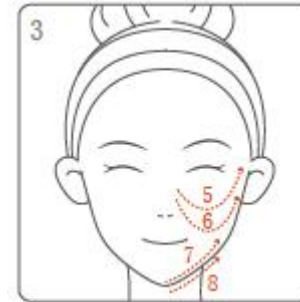

**Right cheek: 50 seconds**

- Apply gel/solution on the right cheek
- Repeat Steps 5-8 (1 set) 3-4 times
- After 50 seconds of use, the device will vibrate briefly, prompting a zone change to the forehead.

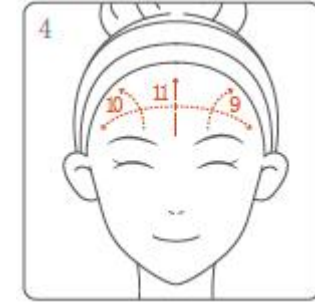

**Forehead: 20 seconds**

- Apply gel/solution on the forehead
- Steps 9-12 make up a set. Repeat 3-4 times.
- After 20 seconds of use, the device will vibrate for an extended period
- Post skincare treatment: Before moving on to the next step in your skincare routine, massage your skin to enhance essence absorption, then wipe off any excess gel or wash it off using a damp cloth.

**Usage Frequency**



Suitable for daily use.

**Warnings**

1. For optimal and consistent energy transmission with minimal loss, it is recommended to utilize Amiro's standard gel or coupling agent.
2. Avoid gliding over the same area more than three times consecutively.
3. Cleanse your face after use.

## After Skincare

---

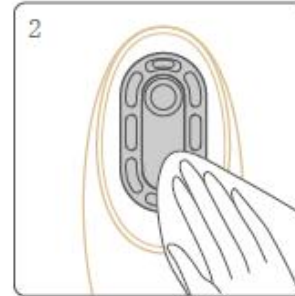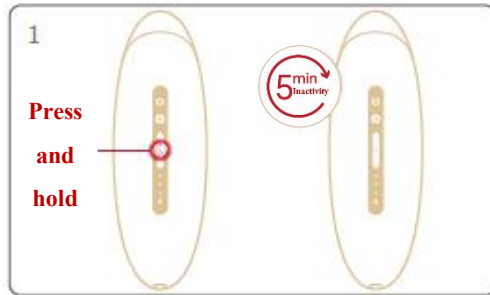

### ① Automatic Shutdown:

- The device will vibrate for an extended duration and automatically switch off after completing skincare in the selected mode.
- The device will automatically switch off after approximately 5 minutes of effective output in Intense Mode.
- The device will automatically switch off after approximately 5 minutes of effective output in Gentle Mode.
- The device will automatically switch off after approximately 2 minutes of effective output in Quick Mode.
- The device will automatically switch off after a certain period of inactivity.

### ② Device Maintenance:

After switching off the device, use a soft, damp cloth or a damp paper towel to wipe the head of the energy box.

#### Caution:

- Avoid cleaning the device with alcohol, alkaline detergents, nail polish remover, or other corrosive liquids.
- Do not clean the film with sharp objects such as nails, toothpicks, needles, etc., and refrain from applying pressure on it.

### ③ Subsequent Skincare:

- Applying a mask for 10-15 minutes after using the device is recommended.
- When the mask treatment is complete, rinse off any remaining essence according to your preference or proceed to the next step in your skincare routine.
- Following the Quick Mode skincare session, make sure to cleanse any leftover gel from your face before proceeding with the next step in your skincare routine.

\*Note: To use the device, build tolerance gradually. Do not use strong skincare products (such as acids/alcohols/whitening products) in the initial stages of using this device to avoid skin discomfort. During treatment, strictly protect your skin from sun exposure and moisturize. Overexposure to the sun can lead to skin sensitivity and photoaging.

## Main Technical Parameters

---

### Conditions for Product Operation

Ambient Temperature: +5°C to +35°C

Relative humidity: 15%-90%

Atmospheric pressure: 70 KPa-106 KPa

Rated input: Type-C 5 V = 3A

Ultrasound Frequency Range: 4.0 MHz ± 10%

### Everyday placement of device




### How to charge the mirror:

Verify that the device is fully charged during the first use

Charging time: ~120 minutes (charging time may vary depending on charging environments)

Battery life: ~60 minutes (battery life may vary depending on usage environments and mode settings)

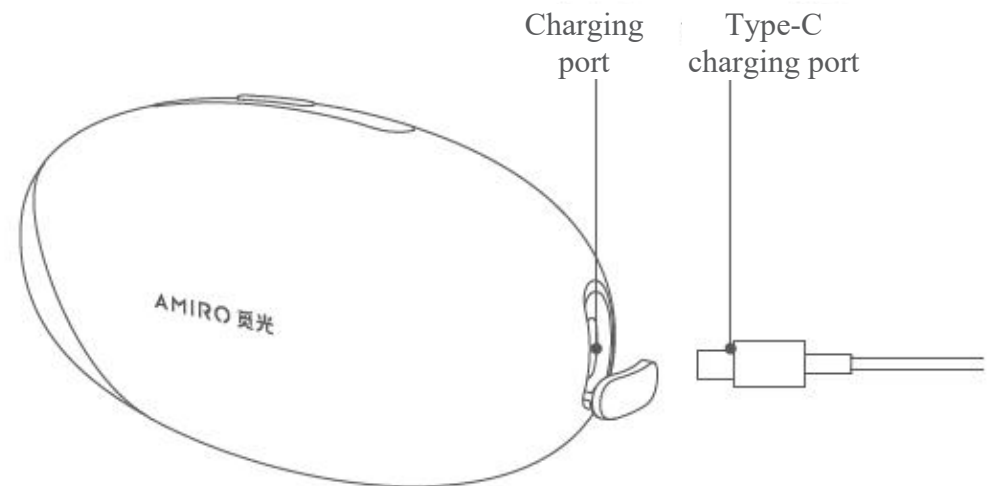



## Function description

<b>Product name</b>		<b>AMIRO Ultrasonic Facial Skin-Tightening Machine</b>
<b>Function pairing</b>	Intense Mode	Ultrasound + Microcurrent + Smart Zone Recognition+ Hydration
	Gentle Mode	Ultrasound + Microcurrent + Smart Zone Recognition+ Skin-Soothing
	Quick Mode	Ultrasound + Microcurrent + Hydration
<b>Materials (Parts contacting human body):</b>	Electrode	17-4 Stainless Steel
	Film	PI Film
	Casing	High Purity Aluminum Alloy
	Base	Silicone + Engineering Plastics

## Product parameters

<b>Product name</b>	AMIRO Ultrasonic Facial Skin-Tightening Machine
<b>Product Model</b>	ABU001
<b>Device Modes:</b>	Intense Mode / Gentle Mode / Express Mode
<b>Device Weight</b>	≈280 g
<b>Ultrasound Frequency</b>	4 MHz ± 10%
<b>Input Voltage</b>	5 V
<b>Rated Input</b>	15 W
<b>Number of Electrodes</b>	8
<b>Microcurrent Frequency</b>	20 KHz±10%
<b>Battery Capacity</b>	3000 mAh
<b>Device Dimensions</b>	135*66*44 mm
<b>Electrode Surface Area</b>	213.07 mm <sup>2</sup>

## Product maintenance

---

- To prevent electrocution, keep the device dry and avoid immersing it or any other parts in water.
- Store the device and its accessories properly after each treatment to prevent damage or loss. Remember to clean the device before putting it away.
- When the device is not in use for an extended period of time, store it in a dry, shaded, and cool place to prevent damage. Avoid contact with hard objects.
- Do not store the device and its accessories in a damp environment for extended periods of time as it can affect the device's performance.
- Before using the device after a long break, make sure to charge it to full, check its functionality, and start at the lowest setting to gradually build tolerance.

## Product storage and transport environment

---

- Ambient Temperature  $-20^{\circ}\text{C} \sim 55^{\circ}\text{C}$
- Ambient Humidity 10%  $\sim$  93%
- Atmospheric Pressure: 70 KPa  $\sim$  106 KPa
- Well-ventilated indoor environment without corrosive gases
- Transportation: Handle with care to prevent impacts, vibrations, and exposure of device to moisture

## Product troubleshooting

FAQ	Possible causes	Solutions
Unable to start the device	Battery depleted	Use the device after it has been charged
	End of battery life	Contact the manufacturer for repair or replacement
	Incorrect operation	Press and hold the power button until the device vibrates for an extended duration and the indicator lights up.
Automatic power off	Not in use for a specific duration	Restart the device
	Device shuts off automatically after powering on	Check the battery indicator light; if it's red, recharge before using again.
	Battery depleted	Charge to full before use

FAQ	Possible causes	Solutions
Device suddenly malfunctions	Smart recognition light flashing red	High temperature, lifespan of energy box expired, unstable connection of energy box
	Not being used for an extended period of time	Device shuts off after vibrating for 5 minutes
	Device malfunction	Please contact the store where you purchased the device or call the customer service hotline at 400-048-0888
Device overheating	Using the device for too long	Let the device rest before using it again.
Heated sensation on the face	Selecting the wrong intensity level	Switch to a gentler level. Start with Level 1 to build tolerance
	Repeated stamping/gliding in the same area	Avoid continuous use on the same area

FAQ	Possible causes	Solutions
Feeling a tingling sensation on the skin	Device not properly fitted against the skin	Ensure the device is in proper contact with the skin
	Choosing the wrong level and mode	Switch to a gentler level. Start with Level 1 to build tolerance Double-click to switch to Gentle Mode for a softer experience Deactivate the microcurrent using Amiro's app
Sensation is too mild	The intensity selected is too low	Switch to higher intensity levels
	Insufficient gel/essence	Increase the gel/essence amount to reduce resistance
	The mask dries out	Apply a fresh one

FAQ	Possible causes	Solutions
Sensation is too mild	Gliding too fast	Decrease the gliding speed and make sure the device glides smoothly and evenly.
	The skin has good tolerance	Skin and nerve sensitivity to electric current varies among individuals. It's a normal phenomenon. Unlock higher energy levels using Amiro's WeChat Mini-Program
Skin discomfort	Redness	Ultrasound is used in this procedure to generate warmth beneath the skin. After the treatment, apply a mask for 10-15 minutes. There will be a temporary recovery period. Seek medical attention promptly if the discomfort persists or does not improve.
	Skin dryness	Apply sunscreen and use hydrating skincare products
Energy box detached	Dropped, caused damage.	Please refer to the instructions on page 32 for installing the energy box. If it fails to operate correctly, kindly reach out to customer service for repair or replacement.










## Main safety properties

---

- 1) The device is rated IPX6 for liquid resistance.
- 2) Input: Type-C 5 V = 3A, 15 W.

## Explanation of graph and symbol

---

Symbol	Explanation	Symbol	Explanation	Symbol	Explanation	Symbol	Explanation
	Warning		Read user manual		Fragile items, handle with care		Keep dry
	Standby/Power		Direct current		Up		Caution
	Environmentally friendly use 5 years						

## Waste disposal declaration

---

### Disposal of used battery:

- Do not attempt to replace the battery on your own during regular use.
- Battery leakage risk: Battery electrolyte is corrosive and should not be in direct contact with the skin or clothing. Promptly wipe the area with a dry cloth if contact occurs. If discomfort persists, seek medical attention immediately.
- Before discarding the product, the battery must be removed from the device:
  - ① Remove the product's casing.
  - ② Disconnect the internal wiring and take out the battery.
- Ensure that the device is powered off and disconnected when removing the battery.
- Safely dispose of the battery.

## Hazardous Substance Label

Component names	Names and concentrations of toxic or hazardous substances or elements in the product					
	Toxic or hazardous substances or elements					
	Lead (Pb)	Mercury (Hg)	Cadmium (Cd)	Hexavalent chromium (Cr(VI))	Polybrominated biphenyls (PBB)	Polybrominated diphenyl ethers (PBDE)
Beauty device PCB	×	○	○	○	○	○
Beauty device casing	○	○	○	○	○	○
Power cable	○	○	○	○	○	○
Base	○	○	○	○	○	○
Electronic components	×	○	○	○	○	○

This table is prepared in accordance with the provisions of SJ/T11364.

○: This indicates that the harmful substance in all homogeneous materials of this component is below the specified limit requirements stated in GB/T-26572-2011.

×: This indicates that the harmful substance in at least one of the homogeneous materials of this component exceeds the specified limit requirements stated in GB/T26572-2011.

The PCB is an internal component of the device that does not come into direct contact with users. The use of such materials will not impact the equipment's safety. This statement serves as a reminder for users to handle the device carefully and adhere to the environmental regulations in their area.



When used as recommended in this manual, this product will last 5 years in an environmentally friendly manner. Within this period, it is safe for use. After this period, it should be properly disposed of and recycled.

## Warranty Policy and Warranty Certificate

---

### Chinese Mainland

- Customers who purchase this product within Chinese Mainland are covered by a 1-year warranty. Warranty periods begin with the date of receipt for online purchases (unless proven otherwise, logistics receipt dates are taken for online purchases; for physical transactions, physical delivery dates are taken for physical transactions).
- Our company delivers products and services in compliance with relevant laws and regulations such as the “E-Commerce Law of the People’s Republic of China”, “Product Quality Law of the People’s Republic of China”, and “Law of the People’s Republic of China on the Protection of Consumer Rights and Interests”. When quality issues arise with in-store purchases within 7 days, contact the seller with a valid receipt or proof of purchase to request a return or exchange. An exchange can be requested from the seller if quality issues arise within 15 days of purchase. For online purchases, you have 7 days from the delivery date to return the product without providing any reason, but it must be in sellable condition. Customers are responsible for paying the shipping charges for returning the product within 7 days, except for returns or exchanges due to quality issues.

### Countries or regions outside Chinese Mainland

- If you purchased this product in a country or region outside the Chinese Mainland, the warranty period, return and exchange policies, and repair procedures will be in accordance with the laws and regulations of that specific country or region. For specific matters, please contact the authorized dealer in your location for resolution.
- Be sure to keep all purchase documentation, including invoices, after purchasing the product. The warranty period will be determined by the production date on the label of the product body if a valid invoice or purchase document is not available during the warranty period.
- Accompany Tech and Shenzhen Accompany Technology Co., Ltd. refer to the same company.



**Please be aware that the following situations are not covered by the free warranty service during the warranty period:**

- Damage caused by improper use, such as incorrect power supply connection, use of inappropriate accessories, failure to follow the instructions, incorrect use, or damage caused by negligence, etc. Repairs or modifications not authorized or approved by our company; other damages caused by force majeure (such as natural disasters, abnormal voltage, etc.).
- Product aging and wear and tear caused by normal use of the product, which does not affect its normal functionality.
- Free servicing does not cover consumable components and add-ons, such as adapters. Delivery charges are not included as well.