



What's Included with Your Collagenius®

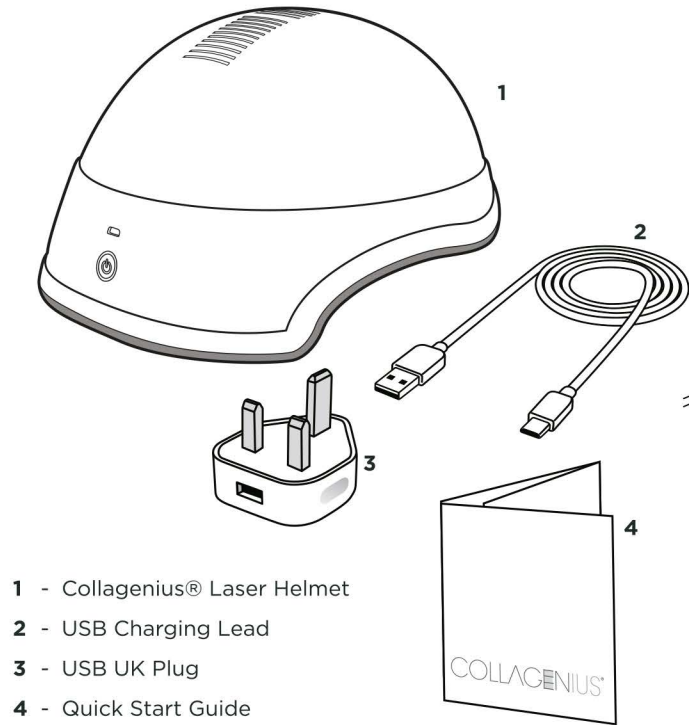

Using Your Collagenius®

Step One

Before your first use your Collagenius® Hair Growth Laser 144 Helmet will need charging. A full charge will take around 4 and a half hours. The power button will illuminate red when charging and go off when complete.

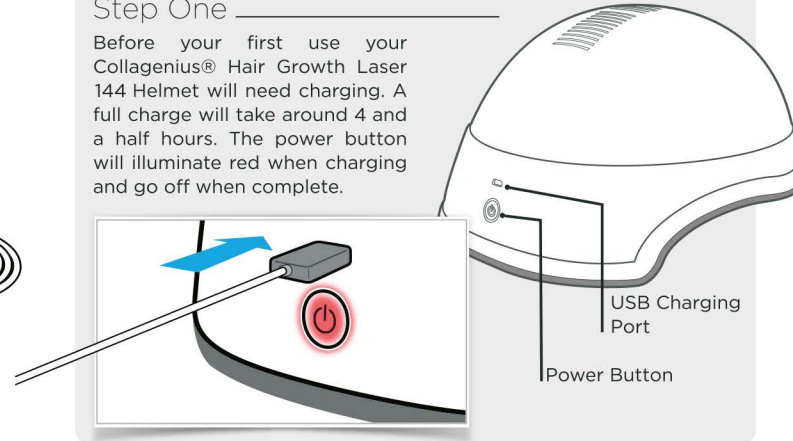


Step Two

Press the power button for 2 seconds, you will hear a 'beep' sound. The power button will turn blue. As a safety feature the laser lights will not illuminate until the helmet is placed on the head.



Step Three

Place the helmet on your head. We recommend wearing it with the power button at the back for comfort. The helmet laser lights will illuminate and start the 25 minute programme, to indicate it is working there is a 'beep' sound every 30 seconds.

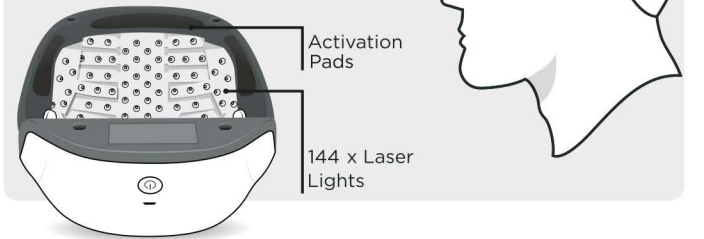

Step Four

If you wish to pause the programme simply remove it from your head. If you remove the helmet during the treatment time the laser lights will switch off for your safety. If you put the helmet back on within 6 minutes the programme will continue as normal. If the helmet is not put back on the head within 6 minutes, there will be a double 'beep' and will turn off automatically ending the programme.

